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DIN EN ISO 4823:2000, Type 1, consistencia Heavy bodied

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PRECISION

DETAIL® HEAVY NEXT

Addition curing silicone with high hardness for pre-impression

INSTRUCTIONS FOR USE

General

Detail® Heavy Next is a precision impression material with high consistency (Heavy Body) for mixing devices, based on Vinyl polysiloxanes (VPS) for addition. Material suggested for double impression technique or double mixing technique combined with Light materials (such as Detail® Light Next Precision line)

Indications

- Impressions on natural teeth for bridges and crowns, inlay and onlay
- Implant impression techniques, if materials with high consistency (Heavy Body) are preferred

Side effects

None known to date. It cannot be excluded that especially predisposed persons might, in singular cases, develop hypersensitivity against single components of the material.

Interactions with other materials

Do not use latex gloves. Avoid contact of the material with organic solvents or liquids containing organic solvents, this could modify its characteristics and cause a lack of precision. Moreover, possible residues of retraction aids can increase the setting time. Do not use the material with condensation curing silicone.

Impression trays and adhesive

Detail® Heavy Next can be used with perforated or non-perforated, metal or resin trays. Ensure using a tray with adequate space available between the tray wall and dental arch. Is it possible to use a tray adhesive for VPS materials in a thin and uniform film, in order to improve the adhesion between Detail® Heavy Next and non-perforated tray. Let the adhesive dry at least 5 minutes on the tray before inserting the impression material.

APPLICATION

1. Charge a sufficient amount of silicone for relining, for exemple Detail® FluidNExt or Detail® Light NExt, in a elastomer syringe and apply in the interested area. At the same time the assistant can put a layer of Detail® Heavy NExt in the prepared tray, paying attention that the mixing tip be always immersed in the material to avoid air bubbles. Loading of the tray and syringing should be completed at the same time, as the higher temperature in the mouth leads to a quicker set of the wash impression material compared to the tray material.
2. Position the tray in the mouth within 1 min. and 30 sec. (Dispenser version) and within 1 min. and 45 sec. (Maxi Dispenser version) from the beginning of the mixing.
3. Remove the impressions after setting of Detail® Heavy NExt, that occurs in not less than 2 min. and 30 sec. (Dispenser version) and in not less than 3 min. and 30 sec. (Maxi Dispenser version). The movement to remove the tray should be parallel to the axis of the prepared teeth.
4. After removal, clean well with water and disinfect the impression. It can be disinfected with standard disinfection solutions. (See "Disinfection and pouring the impression").

Double impression technique

1. Clean and dry the pre-impression done with Detail® Heavy NExt.
2. Then proceed in realizing in the first impression some outlet and remove the imperfections, cleaning and drying carefully.
3. Apply on the pre-impression the silicone for relining (ex. Detail® Fluid NExt or Detail® Light NExt) and position it again in the mouth. Let the material cure.
4. Remove the impressions after setting of Detail® Fluid NExt or Detail® Light NExt that occurs in 3 min. (for Detail® FluidNExt), in 3 min. and 30 sec. (for Detail® Light NExt Normal and Regular), 2 min. (for Detail® Light NExt Fast). The movement to remove the tray should be parallel to the axis of the prepared teeth.
5. After removal, clean well with water and disinfect the impression. It can be disinfected with standard disinfection solutions. (See "Disinfection and pouring the impression"). *Attention: Don't let pass more than 30 min. between the first impression and the second impression.*

Disinfection and pouring the impression

The impression can be disinfected with commercially available immersion solutions. We recommend using a 2% glutaraldehyde solution. Ensure that the disinfectant is indicated for use with addition-cured impression silicones.

Immerse the impression for approx. 10 to 15 minutes. Then rinse under running water for 15 seconds. Keep the impression at normal room temperature (max. 25 °C) and pour the impression no sooner than 30 minutes after removing it from the patient's mouth.

The impression material is compatible with all commercial casts. We recommend using a special purpose hard plaster (ISO 6873, Type 3; e.g. Moldano, Bayer) or an ultra hard dental plaster (ISO

6873, Type 4; e.g. Fujirock II, GC). Impressions can be copperor silver-plated.

ADA recommends following Council on Dental Materials, Instruments and Equipment. Disinfection of impressions. JADA 1991;122(3):110; www.ada.org

PREPARATION OF THE CARTRIDGE

For the version in Dispenser of 50 ml

1. Push up the catch on the back of the dispensing gun and make sure that the plunger be completely inserted into the gun with ratio 2:1.
2. Lift plastic lever on top of the gun and insert the cartridge. Make sure that the notch on the bottom of cartridge aligns with the notch on the dispensing gun.
3. Push down firmly on lever to lock in cartridge.
4. Before attaching a mixing tip, remove the cartridge cap by turning 90° counter clockwise, pull and remove it.
5. When using a cartridge for first time, express a "pea" size amount of material and discard before using a new mixing tip.
6. Line up the notch on the tip with the corresponding notch on the cartridge. Turn it clockwise 90° until it locks into position.
7. The material is dispensed through the mixing tip by pulling the trigger of the dispenser and can be applied directly.
8. Leave a used mixing tip on the cartridge until the next use, it serves as a cap, or remove it and put its cap firmly.

Removal of the Cartridge

Push up catch on the back of the dispensing gun and pull the plunger back. Lift plastic lever on the top of the gun and remove cartridge.


Precautions: To stop the flow, release the trigger. Change the mixing tip before reuse. Remove all trace of hardened paste before replacing the new mixing tip.

Before putting in position the mixing tip, insert the intra-oral tip.

For the Version in Maxi Dispenser of 380 ml

Insert Detail® Heavy NExt cartridge in the mixing device. When doing so observe the user instructions issued by the device manufacturer.

1. Open the locking clip by moving it along the guide to its end position, and then remove the original cap or used mixing tip. (Figure 1-2)

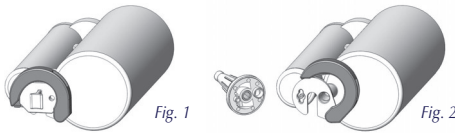2. Now position the new mixing tip and press both openings on the rear side of the mixing tip firmly on to the respective outlet openings in the cartridge. (Figure 3)

3. Once positioned, lock the mixing tip in place with the aid of the locking clip, by moving it along its guide until it surrounds the mixing tip and reaches the end position (Figure 4)

NOTE 1: Always fit a new mixing tip before each mixing procedure! The mixing unit is unable to extrude impression material through set material in the mixing tip.

NOTE 2: Ensure that the mixing tip is fitted securely to the cartridge! Otherwise impression material will be pressed over the mixing tip and may cause contamination.

TECHNICAL DATA

Dispenser Version (Fast Set)

Working time (incl. mixing time) ≤ 1:30 min
Recommended time in mouth ≥ 2:30 min
Strain in compression 1.4 – 2.4 %
Compression set ≤ 0.8 %
Shore-A hardness (after 24 h) ≤ 71

Maxi Dispenser Version (Regular)

Working time (incl. mixing time) ≤ 1:45 min
Recommended time in mouth ≥ 3:30 min
Mixing Ratio (base/catalyst) 5:1
Strain in compression 1.7 – 2.6%
Compression set ≤ 0.7%
Shore-A hardness (after 24 h) ≤ 67

Times refer to normal room temperature (23 °C) and a normal relative humidity (50%). Higher temperatures shorten, lower temperatures increase those times.

Composition

Vinyl polysiloxanes (VPS), silicon dioxide, platinum catalyst, ten-sides and in small amounts stabilizers and additives.

Classification

DIN EN ISO 4823:2000, Type 1, Heavy bodied consistency

Remarks

Store impressions at room temperature and normal humidity.

Storage

Store Detail® Heavy NExt at a room temperature (15° – 25° C). Do not use after the expiry date.

Packaging

Detail® Heavy NExt in Dispenser
2 cartridges of 50 ml each
4 mixing tips

Detail® Heavy NExt in Maxi Dispenser

1 cartridge maxi of 380 ml 10 mixing tips

For dental use only. Keep out of reach of children.

PRECISION products are carefully developed and extensively tested. We therefore guarantee their impeccable quality when used in accordance with their instructions for use. Application, use and processing of the product are beyond our control and are subject to the responsibility of the user.

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DENTALICA S.p.A. – Via Rimini 22 – 20142 Milano – Italia