
mando (6) y aplicar directamente.

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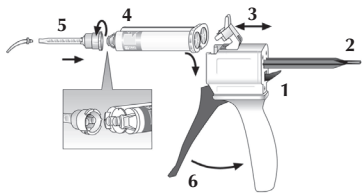2 doppie cartucce da 60 g cad.

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12 beccucci intraorali.

Detail® Light NExt in Dispenser Fast Set

2 doppie cartucce da 60 g cad.

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12 beccucci intraorali.

Detail® Light NExt in Dispenser Regular

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PRECISION

DETAIL® LIGHT NEXT

Addition curing hydrophilic silicone

INSTRUCTIONS FOR USE

General

Detail® Light NExt is a low and medium viscosity thixotropic precision impression material based on addition curing vinyl silicones, supplied in cartridges. It can be applied directly onto the teeth or the tray.

Detail® Light NExt shows hydrophilic properties, thus improving wettability and leading to better reproduction especially of subgingival areas.

Detail® Light NExt is thixotropic. It keeps its position on the tray, but is very fluid under pressure.

For direct extrusion onto teeth intra-oral tips are included.

Indications

Correction impression material for the (two step) double cure the (one step) two-phase impression technique in conjunction with Detail® Putty NExt.

Detail® Light NExt low viscosity comes in handy auto-mixing cartridges in three setting types: Fast Set, Normal Set and Regular.

Composition

Addition curing polysiloxanes, silica, food pigments, additives, catalyst.

Side effects

None known to date.

It cannot be excluded that especially predisposed persons might, in singular cases, develop hypersensitivity against single components of the material.

Interactions with other materials

Do not use with condensation curing silicone impression materials.

Cartridge System

Push up the catch (1) in back of the dispensing gun and completely retract the plunger (2).


Lift plastic lever (3) on top of gun and insert cartridge. It is important to control that notch on bottom of cartridge aligns with the notch on the dispensing gun. Push down firmly on lever (3) to lock in cartridge.

Before attaching a mixing cannula, remove cartridge cap (4) by turning 90° counter clockwise, pull and remove it. Leave an used mixing cannula on the cartridge until it is to be used again, it serves as a cap.

When using a cartridge for first time, express a "pea" size amount of material and discard.

Insert a new mixing cannula (5). Line up the notch on the tip with the corresponding notch on the cartridge. Turn it clockwise 90° until it locks into position.

The material is dispensed through the mixing cannula by pulling the trigger (6) of the dispenser and can be applied directly.

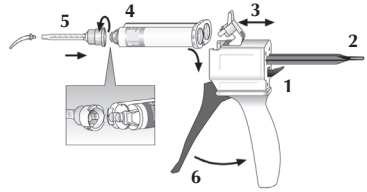Removal of Cartridge

Push up catch (1) in back of the dispensing gun and pull the plunger (2). Lift plastic lever (3) on the top of gun and remove cartridge.

Remarks: To stop the flow, release the trigger.

Change the mixing cannula before re-use.

Remove all trace of hardened paste before replacing the new mixing tip.

Before putting in position the mixing cannula, insert the intra-oral tip.

Impression taking

Detail® Light NExt is loaded directly from the cartridge into the tray or is applied directly onto the teeth. Insert the tray into the mouth no later than 2 minutes 15 seconds (Normal Set) or 1 minute 30 seconds (Fast Set) after application of Detail® Light NExt. Maintain position for at least 3 minutes 30 seconds (Normal Set) or 2 minutes (Fast Set).

NOTE: Times mentioned are valid for 23°C (74°F). Higher temperatures shorten, lower temperatures prolong these times.

Disinfection

Impressions can be disinfected with standard disinfecting solutions.

Casting and plating

The impression should be stored at normal room temperature (max. 25°C). The impression can be poured 30 minutes after removal from the mouth. There are no other time limitations.

Impressions can be copperor silver plated.

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DENTALICA S.p.A. – Via Rimini 22 – 20142 Milano – Italia